



Technomed

Leading the Way

304 MEDICAL GRADE

Model No : 1204 (SLIDING TOP ELECTRIC OT TABLE)



Special Features

This Electromatic sliding top OT Table is suitable for all kind of surgeries, i.e general surgery, vascular surgery, cardio, neurology, urology, gynecology, proctology, laparoscopy, traumatic surgery, plastic surgery and others. The base & column of the table is made of high quality medical stainless steel (Grade-304), non-reflecting surface, antibacterial & easy to clean. We provide with a Top Sliding of 5inch Back/Forward, Height Adjustment, Lateral Tilt, Trendelenburg and Flex/Reflex, Chair Position and Longitudinal Sliding Ttop are realized by Remote Control.

Technical Data

Length of Table	1900 mm
Width of Table Top	520 mm
Height (without Mattress)	Min. 750 to max. 1000 mm
Trendelenburg	+/- 25°
Lateral Tilt	+/- 20°
Back Section	+ /- 60°
Head Section	+/-60°
Leg Section	0°-90°
Floor Locking System	Manual
Flex	+/- 25°
Top Sliding	5" Reverse & 5 " Forward

Salient Features

- Eccentric Pillars
- Stainless Fitting
- Detachable : head, leg & pelvic section
- In built kidney bridge
- Five sectional radio – Translucent Top
- Base & Cover made of attractive, easy to clean Stainless Steel
- Sophisticated mechanics provide smooth, step less articulation of the table top for precise patient positioning
- Standby Battery Back-up with recharging circuit


TABLE POSITIONS

- Height Adjustment – Up & Down
- Lateral Tilt – Left & Right
- Trendelenburg & Reverse Trendelenburg
- Flex & Reflex Position
- Top Sliding

OPTIONAL

- Low Height
- Floor Locking by Remote

STANDARD ACCESSORIES

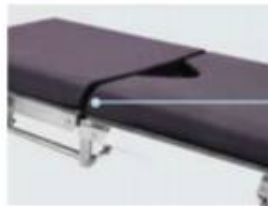

OPTIONAL ACCESSORIES



SPLIT LEG SECTION
OA-1



RESTRAINT BODY STRAP
OA-2



BOLSTER PILLOW
OA-3



INSTRUMENT TRAY
OA-4



FOOT SUPPORT
OA-5



**CERVICAL IN SUPINE /
LAMINECTOMY FRAME**
OA-6



HUMERUS SUPPORT
OA-7



TABLE LENGTH EXTENSION
OA-8



**STEINMANN PIN
CLAMP HOLDER**
OA-9



**SPINE SURGERY USING
CUSHION SPINE FRAME**
OA-10



**ARM TRACTION DEVICE
/ PERINEAL POST**
OA-11



UROLOGY / GYNAE TRAY
OA-12



HORSE SHOE
OA-13



KNEE SUPPORT
OA-14



**FIXATOR FOR
SITTING POSITION**
OA-15



SHOULDER SURGERY
OA-16



**HAND SURGERY
TABLE**
OA-17



WIDTH EXTENSION
OA-18



**ARTHROSCOPY
ATTACHMENT**
OA-19



WRISTLET
OA-20



X-RAY TOP
allows for inserting X-ray film
cassette from head to foot
OA-21



**ORTHOPAEDIC
ATTACHMENT STAND**
OA-22



ACCESSORIES STAND
OA-23

